
Introduction

Thank you for purchasing the Naipo Handheld Percussion Massager with Heat. Please read this user manual carefully and keep it in a safe place for future reference. If you need any assistance, please contact our support team with your product model number and Amazon order number.

Package Contents

- 1 x Percussion Massager
- 3 x Interchangeable massage attachments
- 1 x User Manual

Product Diagram

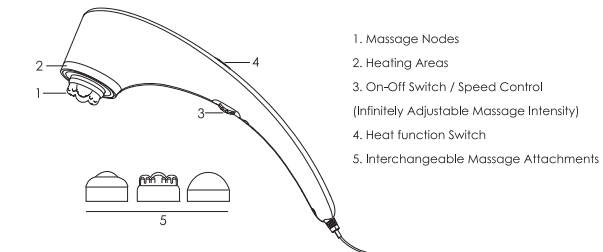Specifications

Dimensions (W*H*D)	95*445*95 mm / 3.74*17.52*3.74 inches
Weight	1000 G
Rated Output	20 W
Rated Voltage	220-240 V ~ 50/60 Hz

Getting Started

- Remove the packaging.
- Check the device, power supply and cables for damage.
- Connect the device to the mains power supply. The device must be switched off.
- Route the cable safely so that there is no risk of anyone tripping over it.

- EN -

Q&A

Problem	Cause	Remedy
Massage head does not move.	Appliance not connected to the power supply.	Connect the plug to the mains supply and switch on the appliance.
	Device does not move.	Switch on the appliance.
	Overheating protection tripped.	Disconnect device from mains power, wait at least 15 minutes, plug in again and switch on.

Product Care

- Disconnect the device before cleaning it.
- Only clean the appliance in the manner specified. Fluid must never be allowed to enter into the device or accessories.
- Do not use the device again until it is completely dry.

Warranty & Customer Support

For questions, assistance, support, or warranty claims, please contact us at the address below that corresponds with your region. Please include your Amazon order number and product model number.

Website: www.naipocare.com

Email: support@naipocare.com

* Please note: Naipo can only provide after sale service for products purchased directly from Naipo. If you have purchased from a different seller, please contact them directly for service or warranty issues.

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Shenzhen Pinlong E-Commerce Co., Ltd.